Repacking Guide



For ComforTrak & FitFlex Table Bed Styles*

*Includes EP-560, EP-860, EP-960, EP-970, NXT-S, 900ia, and 900LX Teeter Inversion Tables.

Please be sure to include all components of the inversion table in your returned order, including the DVD and instructions, otherwise you may not be eligible for a complete refund.

Begin by using the original box and two (2) bags (Large & Small Parts Bag) that came with the inversion table. If you do not have the original box, contact Teeter Customer Service at 800-847-0143. Product returned in an unauthorized packing container may not be eligible for a complete refund.

DO NOT disassemble the product into more components than the parts listed below.

(1) Table Bed

Accessories

6 Roller Hinges

(3) Main Shaft with Ankle Lock System (7) Bolts & Tools Provided

(4) Handles

Start by disconnecting the inversion table into its primary components: Table Bed, A-Frame, Main Shaft with Ankle Lock System, Handles* and Roller Hinges (Figure A). DO NOT disassemble any component into smaller parts.

*For Stretch Max Handle disassembly on the EP-860, EP-960, & EP-970: ONLY untighten the bolt on the spreader arm enough to lift the bottom of the handle off the spreader arm bolt (Figures B & C). DO NOT fully untighten the bolt as this will de-stabilize the A-frame.

Repacking the Table Bed

Remove the accessories (if assembled) from the Table Bed. Place the Table Bed face-down into the box at a slight angle (Figure D).

Repacking the A-Frame

Place the A-Frame into the Large Parts Bag to prevent scratching and then place into the box on top of the Table Bed, with the Crossbar side on the bottom and A-Frame feet at the top of the Table Bed (Figure E). You may need to slightly shift the Table Bed for a proper fit.


Repacking the Ankle Lock System


For T-Pin Ankle Lock:


Close the T-Pin to the last hole setting (Figure F). Then place it Ankle System-side down at an angle next to the Table Bed on top of the plastic bag (Figure C).

(See reverse for EZ-Reach and Deluxe EZ-Reach Ankle Lock System).


NOTE: Figures do not show the Large Parts Bag for easier viewing of box assembly. Inversion Table model below may vary from your own, however same instructions apply.

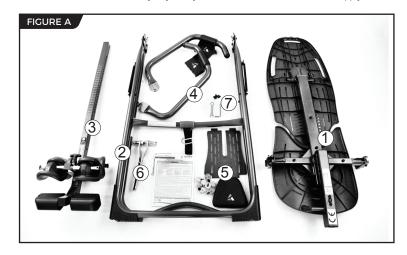

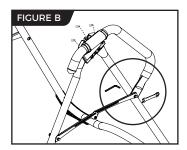

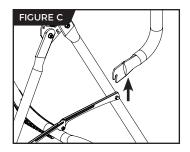



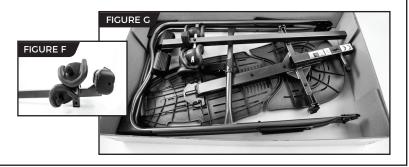




Repacking the Ankle Lock System (Continued)

For EZ-Reach Ankle Lock:

Make sure the EZ-Reach Ankle System is open to the second notch (Figure H). Place Ankle System-side down at an angle next to the Table Bed on top of the plastic bag, with the EZ-Reach Handle tucked under the Crossbar of the A-Frame (Figure I).

For Deluxe EZ-Reach Ankle Lock:

Close the Deluxe Ratchet to the smallest setting (Figure J). Place Ankle System-side down at an angle next to the Table Bed on top of the plastic bag, with the Deluxe EZ-Reach Handle tucked under the Crossbar of the A-Frame (Figure K). You may need to need to slightly lift the A-Frame in order for the Deluxe EZ-Reach handle to slide underneath the crossbar.

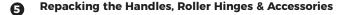

Slide the Handles, Roller Hinges, Accessories, and Product Support materials into the Large Parts Bag (Figure L).


Repacking the Small Parts Bag


Place the Assembly Hardware Kits and Tools into the Small Parts Bag provided. Place Small Parts Bag on top of the Large Parts Bag.

Sealing the Box & Attaching the Label

Seal the large box securely with packing tape, affix the return label and mark your Return Authorization Number (RA #) on the outside panel of the box.

IMPORTANT: Teeter Inversion Tables cannot be returned without a Return Authorization Number. Please contact Teeter Customer Service at (800) 847-014 or email info@teeter.com for assistance.

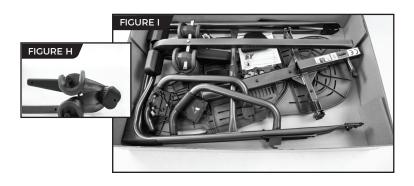

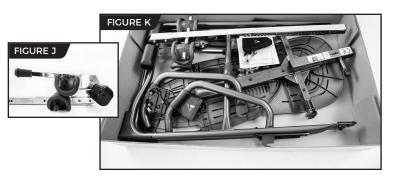





If you have any trouble assembling the equipment, or questions about its use, please contact customer service.

USA: 800-847-0143 or info@teeter.com International: info@teeterintl.com

USA: Teeter 9902 162nd St. Ct. E. Puyallup. WA 98375 Toll Free: 800-847-0143 Fax: 800-847-0188 teeter.com | info@teeter.com International: Teeter International, Ltd.
Gor-Ray House
758 Great Cambridge Rd
Enfield
Middlesex ENI 3GN
United Kingdom
teeterintl.com | info@teeterintl.com